

Product and Specification Guide

Laboratory Case Work and Accessories



Intensa's series of case work offers "Essentials Lab" which can be integrated for all laboratory environments and needs. This flexibility insures adaptations of rapidly changing needs in todays laboratories.

Laboratory Case Work and Accessories



Essentials Lab

Essentials modular, component based system makes installation simple and limits the amount of time needed to turn a space over. Each station can be reconfigured, relocated, or repurposed for ongoing productivity as space and process needs evolve.



Features and Benefits

- * Essentials Lab reconfigurable modular casework can be depreciated in 7 years vs. 39 year depreciation of built in cabinets. Consult tax advisor on specific applications.
- * Essentials Lab casework has a 7 year Limited Warranty. See full warranty details at intensa.net.
- * Essentials Lab casework can be reconfigured and moved at any future time to better accommodate the work that is being accomplished. Built in cabinets cannot be moved or cannot be moved as quickly and easily.
- * Essentials Lab construction provides for units to be shipped fully assembled which means true and accurate alignment of cabinets is assured. They are not shipped in pieces and parts that may mean inaccurate alignment and additional cost to the installers.
- * Base cabinet frames are finished with a durable powder coated finish.
- * All cabinets are finished on the inside and outside. This allows for a finished look as well as prevents any absorption by the substrate. All upper cabinets are designed with white backs for enhanced vision.
- * All cabinets, both base and wall units, are finished on all four sides for sanitary control.
- * All door and drawer edges are finished with heavy duty edge banding. Back edges of the base and wall cabinets are finished to avoid absorption of water or other matter.
- * Essentials Lab uses decorative thermally fused laminate on case sides, back and supporting panels which meets the testing requirements of BIFMA/ANSI laminate requirements.
- * Wall hanging hardware is easily adjustable in the field for proper alignment of cabinets.
- * Intensa, Inc. can provide 2D & 3D drawings plus a quotation at no charge. Customer will need to submit CAD or dimensional drawings.

*****Casework Warning*****

Intensa does not provide fasteners for wall hanging applications. Installers are responsible for furnishing and installing sufficient numbers of the correct fasteners for the field conditions encountered. Installer should always attach wall hanger device into a structural member, or into a heavy duty toggle mechanism, as appropriate for the wall construction and the anticipated load requirements.

Installers are responsible for furnishing and installing adequate fasteners for field conditions encountered.



Laboratory Cabinetry and Case Work



Wall Cabinets



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
1814-SDWC*	18"	14"	30"
2414-DDWC	24"	14"	30"
3014-DDWC	30"	14"	30"
3614-DDWC	36"	14"	30"

*Single door unit please specify Left or Right by denoting "L" or "R"

Wall Cabinets with Tempered Glass



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
1814CA-SDWC*	18"	14"	30"
2414CA-DDWC	24"	14"	30"
3014CA-DDWC	30"	14"	30"
3614CA-DDWC	36"	14"	30"

*Single door unit please specify Left or Right by denoting "L" or "R"

Wall Cabinets Open Front



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
1814-OFWC	18"	14"	30"
2414-OFWC	24"	14"	30"
3014-OFWC	30"	14"	30"
3614-OFWC	36"	14"	30"

Sanitary Slope Top for Wall Cabinets



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
18SLT	18"	14"	4"
24SLT	24"	14"	4"
30SLT	30"	14"	4"
36SLT	36"	14"	4"

Structural Steel Framed Cabinets



Steel Frame

- *Completely welded frame to provide maximum structural stability
- *Durable powder coat finish made to withstand the abuse of an active interior
- *No external fastener connections or finish plugs to allow for tight seamless alignment
- *Threaded inserts for leveling glides or casters
- *All framed units are mobile ready

Finish Panels

- *Produced from dual sided thermally fused laminate available in a wide range of finishes
- *All connections are made internally to allow for easy removal and replacement
- *Exposed edges, doors and drawers are finished with a heavy duty edge banding

Open Front Base Cabinets



ELLA - Includes utility access panel (as shown above)

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYA	18"	24", 30"*	30, 34 and 36"**
ELLA24XXYYA	24"	24", 30"*	30, 34 and 36"**
ELLA30XXYYA	30"	24", 30"*	30, 34 and 36"**
ELLA36XXYYA	36"	24", 30"*	30, 34 and 36"**
ELLA42XXYYA	42"	24", 30"*	30, 34 and 36"**
ELLA48XXYYA	48"	24", 30"*	30, 34 and 36"**
ELLA54XXYYA	54"	24", 30"*	30, 34 and 36"**
ELLA60XXYYA	60"	24", 30"*	30, 34 and 36"**
ELLA66XXYYA	66"	24", 30"*	30, 34 and 36"**
ELLA72XXYYA	72"	24", 30"*	30, 34 and 36"**

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYA	18"	24", 30"*	30, 34 and 36"**
ELLF24XXYYA	24"	24", 30"*	30, 34 and 36"**
ELLF30XXYYA	30"	24", 30"*	30, 34 and 36"**
ELLF36XXYYA	36"	24", 30"*	30, 34 and 36"**
ELLF42XXYYA	42"	24", 30"*	30, 34 and 36"**
ELLF48XXYYA	48"	24", 30"*	30, 34 and 36"**
ELLF54XXYYA	54"	24", 30"*	30, 34 and 36"**
ELLF60XXYYA	60"	24", 30"*	30, 34 and 36"**
ELLF66XXYYA	66"	24", 30"*	30, 34 and 36"**
ELLF72XXYYA	72"	24", 30"*	30, 34 and 36"**

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

Units 48" and wider require vertical divider

Sliding Door Base Cabinets



ELLA - Includes utility access panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA36XXYYB	36"	24", 30"*	30, 34 and 36"**
ELLA42XXYYB	42"	24", 30"*	30, 34 and 36"**
ELLA48XXYYB	48"	24", 30"*	30, 34 and 36"**
ELLA54XXYYB	54"	24", 30"*	30, 34 and 36"**
ELLA60XXYYB	60"	24", 30"*	30, 34 and 36"**
ELLA66XXYYB	66"	24", 30"*	30, 34 and 36"**
ELLA72XXYYB	72"	24", 30"*	30, 34 and 36"**

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF36XXYYB	36"	24", 30"*	30, 34 and 36"**
ELLF42XXYYB	42"	24", 30"*	30, 34 and 36"**
ELLF48XXYYB	48"	24", 30"*	30, 34 and 36"**
ELLF54XXYYB	54"	24", 30"*	30, 34 and 36"**
ELLF60XXYYB	60"	24", 30"*	30, 34 and 36"**
ELLF66XXYYB	66"	24", 30"*	30, 34 and 36"**
ELLF72XXYYB	72"	24", 30"*	30, 34 and 36"**

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

Units 48" and wider require vertical divider

Hinged Door Base Cabinets



ELLA - Includes utility access panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYC***	18"	24", 30"*	30, 34 and 36"***
ELLA24XXYYC	24"	24", 30"*	30, 34 and 36"***
ELLA30XXYYC	30"	24", 30"*	30, 34 and 36"***
ELLA36XXYYC	36"	24", 30"*	30, 34 and 36"***
ELLA42XXYYC	42"	24", 30"*	30, 34 and 36"***
ELLA48XXYYC	48"	24", 30"*	30, 34 and 36"***
ELLA54XXYYC	54"	24", 30"*	30, 34 and 36"***
ELLA60XXYYC	60"	24", 30"*	30, 34 and 36"***
ELLA66XXYYC	66"	24", 30"*	30, 34 and 36"***
ELLA72XXYYC	72"	24", 30"*	30, 34 and 36"***

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYC***	18"	24", 30"*	30, 34 and 36"***
ELLF24XXYYC	24"	24", 30"*	30, 34 and 36"***
ELLF30XXYYC	30"	24", 30"*	30, 34 and 36"***
ELLF36XXYYC	36"	24", 30"*	30, 34 and 36"***
ELLF42XXYYC	42"	24", 30"*	30, 34 and 36"***
ELLF48XXYYC	48"	24", 30"*	30, 34 and 36"***
ELLF54XXYYC	54"	24", 30"*	30, 34 and 36"***
ELLF60XXYYC	60"	24", 30"*	30, 34 and 36"***
ELLF66XXYYC	66"	24", 30"*	30, 34 and 36"***
ELLF72XXYYC	72"	24", 30"*	30, 34 and 36"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

***Single door unit please specify Left or Right by denoting "L" or "R"

Units 48" and wider require double door set

Hinged/Drawer Door Base Cabinets



ELLA - Includes utility access panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA42XXYYD	42"	24", 30"*	30, 34 and 36"**
ELLA48XXYYD	48"	24", 30"*	30, 34 and 36"**
ELLA54XXYYD	54"	24", 30"*	30, 34 and 36"**
ELLA60XXYYD	60"	24", 30"*	30, 34 and 36"**
ELLA66XXYYD	66"	24", 30"*	30, 34 and 36"**
ELLA72XXYYD	72"	24", 30"*	30, 34 and 36"**

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF42XXYYD	42"	24", 30"*	30, 34 and 36"**
ELLF48XXYYD	48"	24", 30"*	30, 34 and 36"**
ELLF54XXYYD	54"	24", 30"*	30, 34 and 36"**
ELLF60XXYYD	60"	24", 30"*	30, 34 and 36"**
ELLF66XXYYD	66"	24", 30"*	30, 34 and 36"**
ELLF72XXYYD	72"	24", 30"*	30, 34 and 36"**

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

Please specify Left or Right by denoting "L" or "R" for drawer placement

Drawer Base Cabinets



ELLA - Includes utility access panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYEZ	18"	24", 30"*	30, 34 and 36"***
ELLA24XXYYEZ	24"	24", 30"*	30, 34 and 36"***
ELLA30XXYYEZ	30"	24", 30"*	30, 34 and 36"***
ELLA36XXYYEZ	36"	24", 30"*	30, 34 and 36"***
ELLA42XXYYEZ	42"	24", 30"*	30, 34 and 36"***
ELLA48XXYYE2Z	48"	24", 30"*	30, 34 and 36"***
ELLA54XXYYE2Z	54"	24", 30"*	30, 34 and 36"***
ELLA60XXYYE2Z	60"	24", 30"*	30, 34 and 36"***
ELLA66XXYYE2Z	66"	24", 30"*	30, 34 and 36"***
ELLA72XXYYE2Z	72"	24", 30"*	30, 34 and 36"***

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYEZ	18"	24", 30"*	30, 34 and 36"***
ELLF24XXYYEZ	24"	24", 30"*	30, 34 and 36"***
ELLF30XXYYEZ	30"	24", 30"*	30, 34 and 36"***
ELLF36XXYYEZ	36"	24", 30"*	30, 34 and 36"***
ELLF42XXYYEZ	42"	24", 30"*	30, 34 and 36"***
ELLF48XXYYE2Z	48"	24", 30"*	30, 34 and 36"***
ELLF54XXYYE2Z	54"	24", 30"*	30, 34 and 36"***
ELLF60XXYYE2Z	60"	24", 30"*	30, 34 and 36"***
ELLF66XXYYE2Z	66"	24", 30"*	30, 34 and 36"***
ELLF72XXYYE2Z	72"	24", 30"*	30, 34 and 36"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

***Z-Specify 4, 3 or 2 drawers (4 drawers shown above)

Units 48" and wider require two drawer sets

Hinged Door/Drawer Base Cabinets



ELLA - Includes utility access panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYG***	18"	24", 30"*	30, 34 and 36"***
ELLA24XXYYG	24"	24", 30"*	30, 34 and 36"***
ELLA30XXYYG	30"	24", 30"*	30, 34 and 36"***
ELLA36XXYYG	36"	24", 30"*	30, 34 and 36"***
ELLA42XXYYG	42"	24", 30"*	30, 34 and 36"***
ELLA48XXYYG	48"	24", 30"*	30, 34 and 36"***
ELLA54XXYYG	54"	24", 30"*	30, 34 and 36"***
ELLA60XXYYG	60"	24", 30"*	30, 34 and 36"***
ELLA66XXYYG	66"	24", 30"*	30, 34 and 36"***
ELLA72XXYYG	72"	24", 30"*	30, 34 and 36"***

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYG***	18"	24", 30"*	30, 34 and 36"***
ELLF24XXYYG	24"	24", 30"*	30, 34 and 36"***
ELLF30XXYYG	30"	24", 30"*	30, 34 and 36"***
ELLF36XXYYG	36"	24", 30"*	30, 34 and 36"***
ELLF42XXYYG	42"	24", 30"*	30, 34 and 36"***
ELLF48XXYYG	48"	24", 30"*	30, 34 and 36"***
ELLF54XXYYG	54"	24", 30"*	30, 34 and 36"***
ELLF60XXYYG	60"	24", 30"*	30, 34 and 36"***
ELLF66XXYYG	66"	24", 30"*	30, 34 and 36"***
ELLF72XXYYG	72"	24", 30"*	30, 34 and 36"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

***Single door unit please specify Left or Right by denoting "L" or "R"

Units 48" and wider require double door/drawer set

Open Desk Base Cabinets



ELLA - Includes utility access panel (as shown above)

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYH	18"	24", 30"*	30, 34 and 36"**
ELLA24XXYYH	24"	24", 30"*	30, 34 and 36"**
ELLA30XXYYH	30"	24", 30"*	30, 34 and 36"**
ELLA36XXYYH	36"	24", 30"*	30, 34 and 36"**
ELLA42XXYYH	42"	24", 30"*	30, 34 and 36"**
ELLA48XXYYH	48"	24", 30"*	30, 34 and 36"**
ELLA54XXYYH	54"	24", 30"*	30, 34 and 36"**
ELLA60XXYYH	60"	24", 30"*	30, 34 and 36"**
ELLA66XXYYH	66"	24", 30"*	30, 34 and 36"**
ELLA72XXYYH	72"	24", 30"*	30, 34 and 36"**

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYH	18"	24", 30"*	30, 34 and 36"**
ELLF24XXYYH	24"	24", 30"*	30, 34 and 36"**
ELLF30XXYYH	30"	24", 30"*	30, 34 and 36"**
ELLF36XXYYH	36"	24", 30"*	30, 34 and 36"**
ELLF42XXYYH	42"	24", 30"*	30, 34 and 36"**
ELLF48XXYYH	48"	24", 30"*	30, 34 and 36"**
ELLF54XXYYH	54"	24", 30"*	30, 34 and 36"**
ELLF60XXYYH	60"	24", 30"*	30, 34 and 36"**
ELLF66XXYYH	66"	24", 30"*	30, 34 and 36"**
ELLF72XXYYH	72"	24", 30"*	30, 34 and 36"**

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

Units 48" and wider require double access panels

Open Desk Return Base Cabinets



ELLA - Includes utility access panel (as shown above)

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLA18XXYYZ	18"	24", 30"*	30, 34 and 36"***
ELLA24XXYYZ	24"	24", 30"*	30, 34 and 36"***
ELLA30XXYYZ	30"	24", 30"*	30, 34 and 36"***
ELLA36XXYYZ	36"	24", 30"*	30, 34 and 36"***
ELLA42XXYYZ	42"	24", 30"*	30, 34 and 36"***
ELLA48XXYYZ	48"	24", 30"*	30, 34 and 36"***
ELLA54XXYYZ	54"	24", 30"*	30, 34 and 36"***
ELLA60XXYYZ	60"	24", 30"*	30, 34 and 36"***
ELLA66XXYYZ	66"	24", 30"*	30, 34 and 36"***
ELLA72XXYYZ	72"	24", 30"*	30, 34 and 36"***

ELLF - Includes full finished back panel

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYYZ	18"	24", 30"*	30, 34 and 36"***
ELLF24XXYYZ	24"	24", 30"*	30, 34 and 36"***
ELLF30XXYYZ	30"	24", 30"*	30, 34 and 36"***
ELLF36XXYYZ	36"	24", 30"*	30, 34 and 36"***
ELLF42XXYYZ	42"	24", 30"*	30, 34 and 36"***
ELLF48XXYYZ	48"	24", 30"*	30, 34 and 36"***
ELLF54XXYYZ	54"	24", 30"*	30, 34 and 36"***
ELLF60XXYYZ	60"	24", 30"*	30, 34 and 36"***
ELLF66XXYYZ	66"	24", 30"*	30, 34 and 36"***
ELLF72XXYYZ	72"	24", 30"*	30, 34 and 36"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

***Z-Specify Left or Right by denoting "L" or "R"

Units 48" and wider require double access panels

Sink Base Cabinets



ELLF - Includes cut out back panel for plumbing

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELLF18XXYY	18"	24", 30"*	30, 34 and 36"***
ELLF24XXYY	24"	24", 30"*	30, 34 and 36"***
ELLF30XXYY	30"	24", 30"*	30, 34 and 36"***
ELLF36XXYY	36"	24", 30"*	30, 34 and 36"***
ELLF42XXYY	42"	24", 30"*	30, 34 and 36"***
ELLF48XXYY	48"	24", 30"*	30, 34 and 36"***
ELLF54XXYY	54"	24", 30"*	30, 34 and 36"***
ELLF60XXYY	60"	24", 30"*	30, 34 and 36"***
ELLF66XXYY	66"	24", 30"*	30, 34 and 36"***
ELLF72XXYY	72"	24", 30"*	30, 34 and 36"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 30", 34" or 36" high

Units 54" and wider require storage doors

Storage Cabinet



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
18XXYY-SDTC***	18"	17.5", 23.5"*	72" and 84"***
24XXYY-DDTC	24"	17.5", 23.5"*	72" and 84"***
30XXYY-DDTC	30"	17.5", 23.5"*	72" and 84"***
36XXYY-DDTC	36"	17.5", 23.5"*	72" and 84"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 72" or 84" high

***Single door unit please specify Left or Right by denoting "L" or "R"

Glass Door Storage Cabinet



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
18XXYYCA-SDTC***	18"	17.5", 23.5"*	72" and 84"***
24XXYYCA-DDTC	24"	17.5", 23.5"*	72" and 84"***
30XXYYCA-DDTC	30"	17.5", 23.5"*	72" and 84"***
36XXYYCA-DDTC	36"	17.5", 23.5"*	72" and 84"***

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 72" or 84" high

***Single door unit please specify Left or Right by denoting "L" or "R"

Open Front Storage



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
18XXYY-OFTC	18"	17.5", 23.5"*	72" and 84"**
24XXYY-OFTC	24"	17.5", 23.5"*	72" and 84"**
30XXYY-OFTC	30"	17.5", 23.5"*	72" and 84"**
36XXYY-OFTC	36"	17.5", 23.5"*	72" and 84"**

* XX-Specify desired finished depth 18" or 24" deep

**YY-Specify desired finished height 72" or 84" high

Wire Shelving



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
XX307C	30"	18", 24"*	74"
XX367C	36"	18", 24"*	74"
XX427C	42"	18", 24"*	74"
XX487C	48"	18", 24"*	74"

* XX-Specify desired finished depth 18" or 24" deep

Open Front Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-OFC	18"	21"	27" or 33"**
SM2421XX-OFC	24"	21"	27" or 33"**
SM3021XX-OFC	30"	21"	27" or 33"**
SM3621XX-OFC	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high


Open Front Mobile Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-OFCM	18"	21"	27" or 33"**
SM2421XX-OFCM	24"	21"	27" or 33"**
SM3021XX-OFCM	30"	21"	27" or 33"**
SM3621XX-OFCM	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

Hinged Door Storage Module

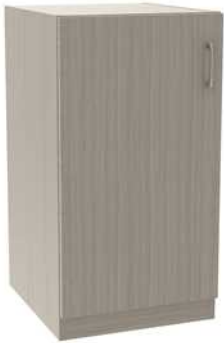

<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-SDC*	18"	21"	27" or 33"**
SM2421XX-SDC*	24"	21"	27" or 33"**
SM3021XX-DDC	30"	21"	27" or 33"**
SM3621XX-DDC	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

*Single door unit please specify Left or Right by denoting "L" or "R"

Hinged Door Mobile Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-SDCM*	18"	21"	27" or 33"**
SM2421XX-SDCM*	24"	21"	27" or 33"**
SM3021XX-DDCM	30"	21"	27" or 33"**
SM3621XX-DDCM	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

*Single door unit please specify Left or Right by denoting "L" or "R"

Hinged Door/Drawer Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-SDRW*	18"	21"	27" or 33"**
SM2421XX-SDRW*	24"	21"	27" or 33"**
SM3021XX-DDRW	30"	21"	27" or 33"**
SM3621XX-DDRW	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

*Single door unit please specify Left or Right by denoting "L" or "R"

Hinged Door/Drawer Mobile Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-SDRWM*	18"	21"	27" or 33"**
SM2421XX-SDRWM*	24"	21"	27" or 33"**
SM3021XX-DDRWM	30"	21"	27" or 33"**
SM3621XX-DDRWM	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

*Single door unit please specify Left or Right by denoting "L" or "R"

Drawer Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-EZ	18"	21"	27" or 33"**
SM2421XX-EZ	24"	21"	27" or 33"**
SM3021XX-EZ	30"	21"	27" or 33"**
SM3621XX-EZ	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

***Z-Specify 4, 3 or 2 drawers (4 drawers shown above)

Drawer Mobile Storage Module



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
SM1821XX-EZM	18"	21"	27" or 33"**
SM2421XX-EZM	24"	21"	27" or 33"**
SM3021XX-EZM	30"	21"	27" or 33"**
SM3621XX-EZM	36"	21"	27" or 33"**

**XX-Specify desired finished height 27" or 33" high

***Z-Specify 4, 3 or 2 drawers (4 drawers shown above)

Adjustable Table Frame



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELT2424A	24"	23.5"	26-36"
ELT3024A	30"	23.5"	26-36"
ELT3624A	36"	23.5"	26-36"
ELT4224A	42"	23.5"	26-36"
ELT4824A	48"	23.5"	26-36"
ELT5424A	54"	23.5"	26-36"
ELT6024A	60"	23.5"	26-36"
ELT6624A	66"	23.5"	26-36"
ELT7224A	72"	23.5"	26-36"
ELT2430A	24"	29.5"	26-36"
ELT3030A	30"	29.5"	26-36"
ELT3630A	36"	29.5"	26-36"
ELT4230A	42"	29.5"	26-36"
ELT4830A	48"	29.5"	26-36"
ELT5430A	54"	29.5"	26-36"
ELT6030A	60"	29.5"	26-36"
ELT6630A	66"	29.5"	26-36"
ELT7230A	72"	29.5"	26-36"

Adjustable Table Frame with Shelf



<u>Item</u>	<u>Width</u>	<u>Depth</u>	<u>Height</u>
ELT2424B	24"	23.5"	26-36"
ELT3024B	30"	23.5"	26-36"
ELT3624B	36"	23.5"	26-36"
ELT4224B	42"	23.5"	26-36"
ELT4824B	48"	23.5"	26-36"
ELT5424B	54"	23.5"	26-36"
ELT6024B	60"	23.5"	26-36"
ELT6624B	66"	23.5"	26-36"
ELT7224B	72"	23.5"	26-36"
ELT2430B	24"	29.5"	26-36"
ELT3030B	30"	29.5"	26-36"
ELT3630B	36"	29.5"	26-36"
ELT4230B	42"	29.5"	26-36"
ELT4830B	48"	29.5"	26-36"
ELT5430B	54"	29.5"	26-36"
ELT6030B	60"	29.5"	26-36"
ELT6630B	66"	29.5"	26-36"
ELT7230B	72"	29.5"	26-36"

Essentials Lab



Workstations and Accessories



Concept Workstation Frame



	<u>Manual</u>	<u>Handcrank</u>	<u>Motor</u>
30x30	E827281-49	E824551-49	E825361-49
30x36	E827283-49	E819146-49	E825156-49
24x40	10049015P	10149015	10549015
30x40	10049030P	10149030	10549030
36x40	10049062P	14-1014901	10549062P
24x48	10049016P	10149016	10549016
30x48	US-10049031P	10149031	10549031
36x48	10049063P	10149063	10549063P
24x60	10049017P	10149017	10549017
30x60	US-10049032P	10149032	10549032
36x60	10049064P	10149064	10549064P
24x72	10049018P	10149018	10549018
30x72	US-10049033P	10149033	10549033
36x72	10049065P	10149065	10549065P
24x80	10049019P		10549019
30x80	10049034P		10549034
36x80	10049066P		10549066P
24x96	14-10049436		
30x96	E827370-49		14-10349249
36x96	14-1004968		14-10349283

Manual Frame - adjustable height 26.62-44.34 inches

Handcrank Frame - adjustable height 26.82-40.60 inches

Motor Frame - adjustable height 27.80-43.55 inches

*Standard powder coat finish light gray (ESD Powder Coat)

*Tops ordered separately

Cornerstone Workstation Frame



Frame

CS712-49	Cornerstone frame, 30" x 48" 27.56" x 47.24" x 24.72-37.64"
CS715-49	Cornerstone frame, 30" x 60" 27.56" x 59.06" x 24.72-37.64"
CS718-49	Cornerstone frame, 30" x 72" 27.56" x 70.87" x 24.72-37.64"
E903146-49	Cornerstone frame, 30" x 80" x 27.56" x 94.49" x 24.72-37.64"
E827221-49	Cornerstone frame, 30" x 96" 27.56" x 94.49" x 24.72-37.64"

Depth Kit

CSD60-49	Cornerstone depth kit 24"*
CSD90-49	Cornerstone depth kit 36"*

*Cornerstone depth kits should be ordered with a 30" deep frame for 24" and 36"d tables

Lower Shelf

LSSCS120-49	Lower steel shelf 40.94"x18.1" Cornerstone
LSSCS150-49	Lower steel shelf 53.15"x18.1" Cornerstone
LSSCS180-49	Lower steel shelf 64.96"x18.1" Cornerstone

*Standard powder coat finish light gray (ESD Powder Coat)

*Tops ordered separately

Caster sets for Concept and Cornerstone frames

834734-00	Caster set 3.93" NonESD
14-7000025	Roll & Set caster set



Workstation Accessories

Uprights (53"H)

1 x M30 Upright module	90049001P
1 x M36 Upright module	90149001P
1 x M48 Upright module	US-9074901P
1 x M60 Upright module	US-9074908P
1 x M72 Upright module	US-9074909P
2 x M30 Upright module	US-90349001P
2 x M36 Upright module	US-90449001P



Tiltable Steel Shelf (includes one lip that can be mounted on front or back)

Shelf M30	28.34" x 11.81"	TTSM7503-49
Shelf M30	28.34" x 15.74"	TTSM7504-49
Shelf M36	34.25" x 11.81"	TTSM9003-49
Shelf M36	34.25" x 15.74"	TTSM9004-49
Shelf M48	45.94" x 11.81"	TTSM1203-49
Shelf M48	45.94" x 15.74"	TTSM1004-49
Shelf M60	57.76" x 11.81"	TTSM15003-49
Shelf M60	57.76" x 15.74"	TTSM15004-49
Shelf M72	69.80" x 11.81"	TTSM18003-49
Shelf M72	69.80" x 15.74"	TTSM18004-49



Lighting

Light/Balancer rail M48	US-91849010P
Light/Balancer rail M60	US-91849006P
Light/Balancer rail M72	US-91849008P
36" LED light, balancer rail hardware	14-95035172
36" LED light, magnetic attachment	14-95035175
48" LED light, balancer rail hardware	14-95035173
48" LED light, magnetic attachment	14-95035176



Bin Rails

Bin rail M30	29.01"	US-859151-49
Bin rail M36	35.43"	US-859152-49
Bin rail M48	47.04"	US-859157-49
Bin rail M60	58.89"	US-859154-49
Bin rail M72	70.66"	US-859155-49



Workstation Accessories

Power Rails

M30 horizontal power rail	14-9114905
M36 horizontal power rail	14-9114910
M48 horizontal power rail	14-9114915
M60 horizontal power rail	14-9114919
M72 horizontal power rail	14-9114923
Vertical power unit with clamp	14-9114928



IT/Tech Support

CPU holder, leg mount	860552-49
CPU holder, motor leg mount	860740-49
LCD Holder with double arm	US-93049002P
Keyboard Tray for LCD holder	860941-49
Under surface keyboard tray	92549001P
Height adjustable LCD arm	14-94049440
CPU holder for Cornerstone	860553-49



Wall Mount Accessories

Oval rails and Upright tubes

Oval horizontal rail 96"	890580-49
Oval horizontal rail 48"	E816366-49
Upright tube for oval rail 47.24"	100071-49
Upright tube for oval rail 78.74"	100070-49
End caps and plug set	14-1094902
Extension for oval rail	110913-51
End caps and plug set	010014-00



Preconfigured Reagent Rack



Top Mounted Reagent rack

	30"W	36"W	48"W	60"W	72"W
Shelf Bay 31.10"	14-9564907	14-9564908	14-9564909	14-9564938	14-9564939
Extension 31.10"	14-9564910	14-9564911	14-9564912	14-9564940	14-9564941
Shelf Bay 39.37"	14-9564942	14-9564943	14-9564944	14-9564945	14-9564946
Extension 39.37"	14-9564947	14-9564948	14-9564949	14-9564950	14-9564951
Shelf Bay 47.24"	14-9564952	14-9564953	14-9564954	14-9564955	14-9564956
Extension 47.24"	14-9564957	14-9564958	14-9564959	14-9564960	14-9564961


Shelf Bay Includes:

- 2 pcs Uprights w/brackets
- 1 pc Top cross rail
- 2 pcs Corners

Shelf Bay Extension Includes:

- 1pc Upright w/bracket
- 1 pc Top cross rail
- 1 pc T bracket

Utility Support Upright System

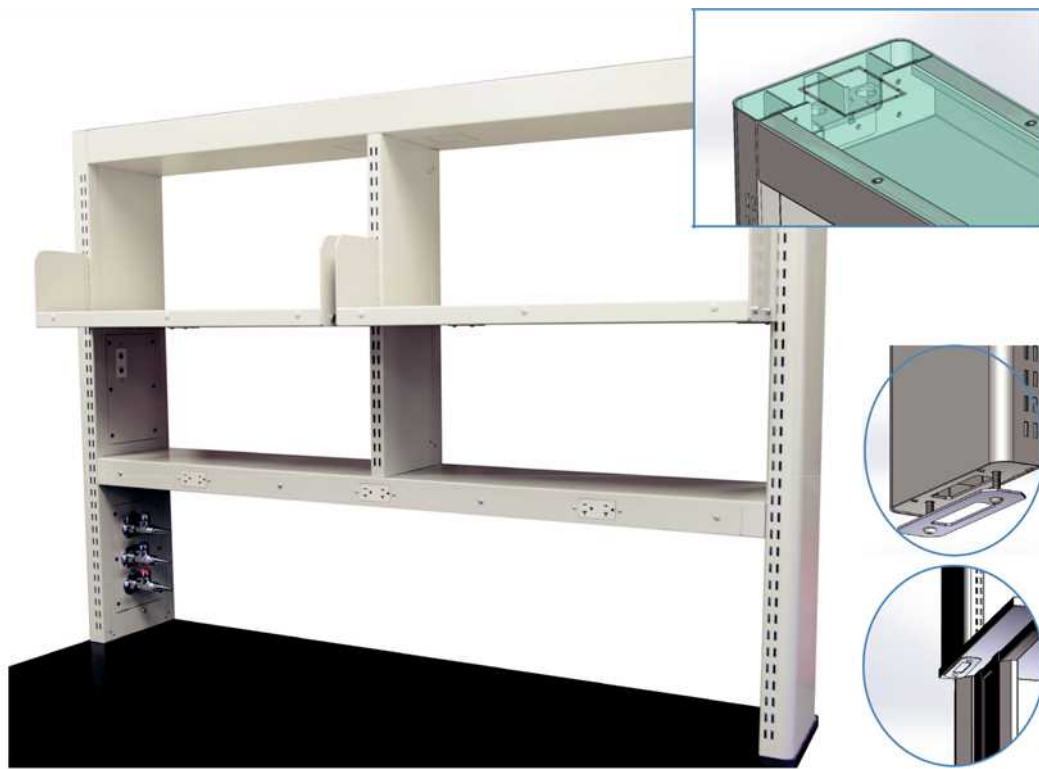

USU Single Bay Kits

M48 Single Bay Counter Mounted USU

USUM1200-49

M60 Single Bay Counter Mounted USU

USUM1500-49

M72 Single Bay Counter Mounted USU

USUM1800-49

Single Bay Starter Kits Include (2) vertical USU uprights, blank service panels with attachment hardware*, fixed horizontal utility channel/shelf with blank front panel*, and upper horizontal cable channel with cover.

*Blank panels do not include cut-outs for connecting services. To add services, please select the panel according to the services requested.

Shelves for Single Bay Kits

M48 x 12"D USU Shelf with Bookend Brackets

USUSHM1203-49

M48 x 16"D USU Shelf with Bookend Brackets

USUSHM1204-49

M60 x 12"D USU Shelf with Bookend Brackets

USUSHM1503-49

M60 x 16"D USU Shelf with Bookend Brackets

USUSHM1504-49

M72 x 12"D USU Shelf with Bookend Brackets

USUSHM1803-49

M72 x 16"D USU Shelf with Bookend Brackets

USUSHM1804-49

*Shelves do not contain utility support. Utility support is housed in vertical USU uprights and horizontal channel only.

Utility Support Upright System

Horizontal Utility Channels

M48 Horizontal utility channel	USUHCM1200-49
M60 Horizontal utility channel	USUHCM1500-49
M72 Horizontal utility channel	USUHCM1800-49

Horizontal utility channel includes a blank front panel and bottom cover. Additional panels available for front-facing utility support. Horizontal utility channel mounts to single bay USU system. Adding a middle divider to this channel creates a double bay system.

Service Panels with Cut-Outs for Horizontal Utility Channels

M48 Horizontal Panel with 3 cut-outs for Decora Duplex	USUFPM1203-49
M48 Horizontal Panel with 4 cut-outs for Decora Duplex	USUFPM1204-49
M48 Horizontal Panel with 5 cut-outs for Decora Duplex	USUFPM1205-49
M60 Horizontal Panel with 3 cut-outs for Decora Duplex	USUFPM1503-49
M60 Horizontal Panel with 4 cut-outs for Decora Duplex	USUFPM1504-49
M60 Horizontal Panel with 5 cut-outs for Decora Duplex	USUFPM1505-49
M72 Horizontal Panel with 3 cut-outs for Decora Duplex	USUFPM1803-49
M72 Horizontal Panel with 4 cut-outs for Decora Duplex	USUFPM1804-49
M72 Horizontal Panel with 5 cut-outs for Decora Duplex	USUFPM1805-49

*Horizontal service panels with cut-outs replace blank horizontal channel panels in starter kits.

Middle Divider and Shelves for USU Double Bay

*Middle Vertical Divider (slotted)	USUDIV-49
M30 x 12"D USU Shelf with Bookend Brackets	USUSHM7503-49
M30 x 16"D USU Shelf with Bookend Brackets	USUSHM7504-49
M36 x 12"D USU Shelf with Bookend Brackets	USUSHM9003-49
M36 x 16"D USU Shelf with Bookend Brackets	USUSHM9004-49

*Middle vertical divider is required for attachment of M30 and M36 shelves.

**Divider and shelves do not contain utility support. Utility support is housed in vertical USU uprights and horizontal channel only.

Service Panels with Cut-Outs for Vertical USU Uprights

Service Panel with (3) 1" Holes for Services	USUSP11-49
Service Panel with 1 Duplex Power Cut-Out	USUSP12-49
Service Panel with 2 Duplex Power Cut-Outs	USUSP13-49
Service Panel with 2 Duplex Power Cut-Outs and 1 Data Cut-Out	USUSP14-49

Essentials Lab



More Space Post, Rails and Accessories



More Space Post, Rails and Accessories



Wall Mount



Bolt Mount



Clamp Mount



Grommet Mount

Wall Mount

961+0509

17-1/2"H Wall Mount Post

961+0519

33-1/4"H Wall Mount Post

Bolt Mounted

961+0109

17-1/2"H Bolt Mount Post

961+0139

21-1/2"H Bolt Mount Post

961+0119

33-1/4"H Bolt Mount Post

Clamp Mount

961+0609

17-1/2"H Clamp Mount Post

961+0639

21-1/2"H Clamp Mount Post

Grommet Mount

961+0409

17-1/2"H Grommet Mount Post

961+0439

21-1/2"H Grommet Mount Post

961+0419

33-1/4"H Grommet Mount Post

*Post are 2" diameter extruded aluminum

*Includes hardware and wire management clips

More Space Post, Rails and Accessories



SlatWall-Element

- 973+0805 Slatwall element 800mm (31.5") wide
- 973+1005 Slatwall element 1000mm (39.4") wide
- 973+1205 Slatwall element 1200mm (47.25") wide
- 973+1605 Slatwall element 1600mm (63") wide



Pencil cup



Tray



Memo box



Multi-Purpose box



Multi-storage board

Rail Accessories

- 725+0019 Tray B4 portrait format Penda
- 725+0029 Tray B4 landscape format Penda
- 725+0039 Memo box Penda
- 725+0049 Small multi-purpose box Penda
- 725+0059 Large multi-purpose box Penda
- 725+0069 Pencil cup Penda
- 725+0079 Phone holder Penda
- 725+0089 Tablet holder Penda
- 725+0099 Multi storage board Penda
- 725+0129 File holder Penda



Phone/Tablet holder



File holder

More Space Post, Rails and Accessories



Monitor Support Arms

962+0119	TSS folding arm I
963+0119	TSS support carriage
964+0119	TSS folding arm II
965+0119	TSS folding arm III
967+0119	TSS keyboard holder



LiftTEC Arms

980+1089	LiftTEC Arm I - No forearm Column mount 6-18lbs
980+1159	LiftTEC Arm I - No forearm Column mount 15-33lbs
980+2089	LiftTEC Arm II - Angled forearm Column mount 6-8lbs
980+2159	LiftTEC Arm II - Angled forearm Column mount 15-33lbs
980+3089	LiftTEC Arm III - Straight forearm Column mount 6-8lbs
980+3159	LiftTEC Arm III - Straight forearm Column mount 15-33lbs
981+1089	LiftTEC Arm I - No forearm SlatWall mount 6-18lbs
981+2089	LiftTEC Arm II - Angled forearm SlatWall mount 6-8lbs
981+3089	LiftTEC Arm III - Straight forearm SlatWall mount 6-8lbs
931+1089	LiftTEC Arm I - No forearm Wall mount 6-18lbs
931+2089	LiftTEC Arm II - Angled forearm Wall mount 6-8lbs
931+3089	LiftTEC Arm III - Straight forearm Wall mount 6-8lbs

Sinks, Faucets and Accessories



Stainless drop in sinks

SL1515-B-GR	Single bowl with ledge 15"x15"x6"
SL1921-B-GR	Single bowl with ledge 19"x21"x7"
DL1925-B-GR	Double bowl with ledge 19"x25"x7"
DL2133-B-GR	Double bowl with ledge 21"x33"x7"



Epoxy drop in sinks

CRD10	Drop in sink 16"x8"x7"
CRD15	Drop in sink 16"x12"x8"
CRD25ADA	Drop in sink 18"x15"x5"
CRD25	Drop in sink 18"x15"x8"
CRD30	Drop in sink 18"x15"x11"
CRD55ADA	Drop in sink 25"x15"x5"
CRD55	Drop in sink 25"x15"x10"
CRD59	Drop in sink 28"x15"x12"
CRC36	Epoxy 3"x6"
CRC39	Epoxy 3"x9"

Sinks, Faucets and Accessories



Mixing Faucets

L2224

6" Swing gooseneck with paddle handles aerator nozzle

L2224VB

6" Swing gooseneck with paddle handles with vacuum breaker

L412

6" Swing gooseneck with 4-arm handles serrated nozzle

L412VB

6" Swing gooseneck with 4-arm handles with vacuum breaker

Single Faucets

L611

6" Swing gooseneck with 4-arm handle serrated nozzle

L611VB

6" Swing gooseneck with 4-arm handle with vacuum breaker

Pure Water Faucets

L681

6" Swing gooseneck with 4-arm handle serrated nozzle

L681VB

6" Swing gooseneck with 4-arm handle with vacuum breaker



EW1022

Eyewash with drench hose
Deck mounted

Sinks, Faucets and Accessories



Stainless Peg Racks

- V3030 30"x30" wall mounted with 50 pegs
- V2430 24"x30" wall mounted with 32 pegs



Needle Valve Deck Mounted

- L2880-131WSA Deck mounted single valve
- L2880-132AWSA Deck mounted double valve 90 degrees
- L2880-132SWSA Deck mounted double valve 180 degrees

Needle Valve Panel Mounted

- L4880F-225WSA Panel mounted single valve
- L2880-141WSA Panel mounted double valve

Accessories

<u>Item</u>	<u>Item Number</u>	<u>Description</u>
Stainless Grommet	SGR-04	4" Flush Mount Grommet
	SGR-06	6" Flush Mount Grommet
	SGR-08	8" Flush Mount Grommet
Wire Grommet	GR01	1-5/8" Black Grommet with Cover
	GR03	3" Black Grommet with Cover
Task Light	TL01	Cool White LED Task Light
Interior Shelf	IS-XXYY	Adjustable Interior Shelf
Vertical Partition	VP-XXYY	Cabinet Vertical Partition
Work Surface Bracket	WSB-XX	Wall Bracket for Counter
Filler Panel	FP-XXYY	Filler Panel and Hardware
Casters	FC-01	HD Dual Wheel Furniture Casters (Full Swivel-2 Brakes) Black
		2 Outlet flip up power unit (Ground Plug)
In Surface Power	BE02137	2 Outlet flip up power unit (Ground Plug)
	BE080040	2 Outlet flip up power unit (Hardwire)

Locks

<u>Item</u>	<u>Item Number</u>	<u>Description</u>
Hinged Door Lock	HL	Key Lock with Changeable Core
Drawer Lock	SDL	Key Lock with Changeable Core
Gang Lock	GL	Ganged Key Lock for Drawer Stack (Changeable Core)
Combo Lock	CL	3 Digit Combination Lock

*For Any Other Style Lock Requirements Please Contact Intensa

Standard Laminate Finishes



Charleston Oak
(CHO)



Shaker Cherry
(SHA)



Wild Cherry
(WCH)



Williams Oak
(WMO)



Norway
(NWX)



Ganache
(GAN)



Solar Oak
(SOK)



Talas Cherry
(TCY)



Windsor Mahogany
(WMY)



Aria
(ARA)



Black*
(BLK)



Folkstone Grey*
(FGY)



London Grey*
(LGY)




Icicle*
(ICE)

Standard Powder Coat Finishes



Black
(BLK)

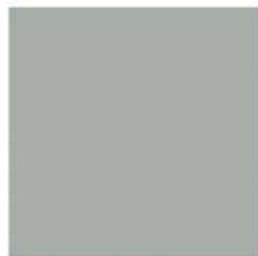

London Grey
(LGY)



Folkstone Grey
(FGY)



Icicle
(ICE)



Silver
(SLR)

Counter Top Options

Phenolic Resin

Phenolic resin is made by fusing layers of craft paper saturated in resin. The resin-soaked layers combine to form one solid piece of countertop that's designed for longevity. Available in 1" and 3/4" thickness. Phenolic resin countertops are highly chemical resistant and extremely durable. These are great countertops for labs with continuous heat exposure up to 350 degrees Fahrenheit. Phenolic resin can be easily modified in the field.

Epoxy Resin

Epoxy resin has long been the ultimate standard countertop material for most laboratories. It's non-porous, making it the perfect choice for labs dealing with caustic chemicals, and doesn't absorb water so is often found around sinks and faucets. Epoxy resin is made by combining a variety of materials—silica, filler material, a hardener, and epoxy resin—and oven-curing the materials in one solid piece. Available in 1" thickness.

Stainless Steel

Obviously one of the sleekest choices, stainless steel countertops are highly durable and made from one material – stainless steel. Many labs opt for stainless steel when they're looking for visual continuity that pairs with overall function. Insanely durable, stainless steel has a lot of benefits, but it might not be the countertop choice for everyone.

Maple Block

Also known as butcher block, wood countertops are made by laminating hard maple planks together to create a two-inch thick top. This is one of the more affordable choices and can be used in just about any application where chemicals aren't a big concern.

High Pressure Laminate

One of the more affordable options, plastic laminate is made up of a combination of phenolic resin, melamine, and wood fibers. It's a great all-around countertop that's best used in dry working conditions like technology labs. High pressure laminate is also available in a chemical resistant version. All high pressure counters are finished with a heavy duty edge treatment.

Intensa, Inc. | P.O. Box 5981 | High Point, NC 27262 | USA
+1.888.636.2660 | sales@intensa.net | intensa.net

Casework Manufactured in the USA
GF2000221RevA20

